

cleanoz[®]

nasal aspirator kit, aspirador nasal
kit mouche bébé

MBO02



QUICK, RÁPIDO, RAPIDE

EFFECTIVE, EFICAZ, EFFICACE

HYGIENIC, HIGIÉNICO, HYGIENIQUE

EASY TO USE, FÁCIL, FACILE A UTILISER

GENTLE, SUAVE, DOUX


SAFE, SEGURO, SUR



ubimed[®]
USER MANUAL

USO MANUAL
NOTICE D'UTILISATION



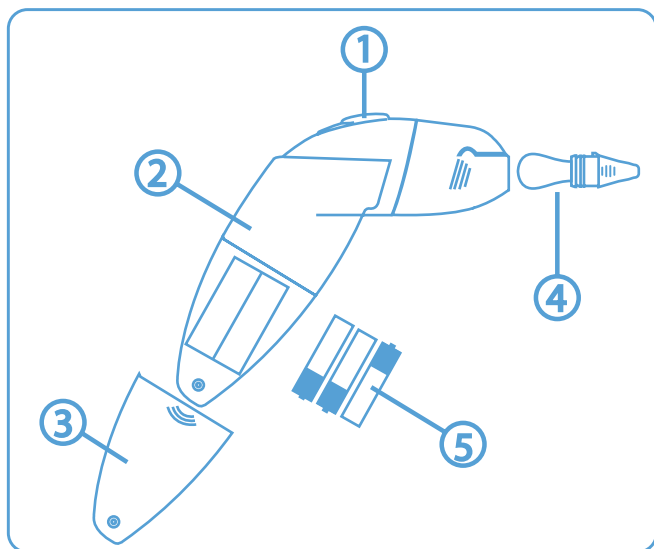

GENERAL DESCRIPTION

- ① Switch On-Off
- ② Main Unit and transparent tank
- ③ Battery door
- ④ Disposable Reservoir Nozzle tip
- ⑤ Batteries

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE



CLEANOZ® has been carefully engineered for both efficacy and safety. However, basic safety precautions should always be followed when using any electrical device, especially intended for baby's care.

WARNING

Use this product only for its intended use as described in this booklet.

Do not use attachments not recommended by Ubimed®.

RISK OF EXPLOSION IF BATTERY IS REPLACED BY AN INCORRECT TYPE.

REPLACE ONLY WITH ALKALINE BATTERIES AAA LR03.


For any question regarding Cleanoz®, please contact Ubimed Customer Service Department online at customer-service@ubimed.com

or our Call Center: 1-866-5-UBIMED, US only (9 AM to 5 PM PST, Monday to Friday).

MEDICAL WARNING

CAUTION: THE DISPOSABLE RESERVOIR NOZZLE PIECE CONTAINS A LATEX BAG WHICH MAY CAUSE ALLERGENIC REACTIONS.

- 1) Wash your hands before and after using Cleanoz®.
- 2) Do not push nozzle deep into baby's nostril. Use Cleanoz® in the recommended direction.

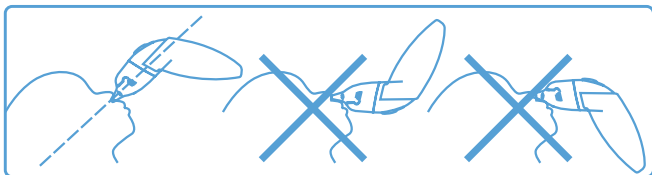


- 3) For the baby's security, aspiration time is limited to few seconds.
- 4) Do not use Cleanoz® in case of nasal irritation or hemorrhage.
- 5) In the case of infection, fever, nasal irritation or hemorrhage, please consult a doctor.
- 6) Do not use Cleanoz® without disposable nozzle piece or with a deteriorated bag.
- 7) Use Cleanoz® only with Ubimed's disposable nozzle piece.
- 8) Change the disposable nose piece after each use.

SAVE THESE INSTRUCTIONS


INSTRUCTIONS FOR USE

Cleanoz® is an electric nose cleaner for babies


1) For the first use: Open the battery compartment by removing the cover from the nose cleaner (use a screwdriver to unscrew). Insert three AAA LR3 Alkaline batteries then shut the compartment: Use only Alkaline batteries and respect the polarity.

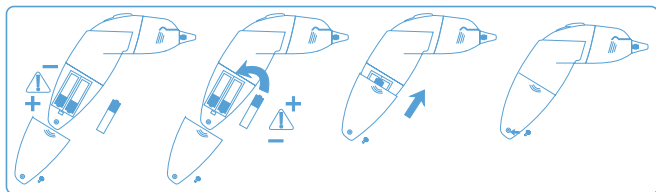

2) Insert a new disposable reservoir nozzle piece by aligning the upper silicone guide nodule with the excised edge of the transparent tank. Please be careful not to damage the bag upon insertion.

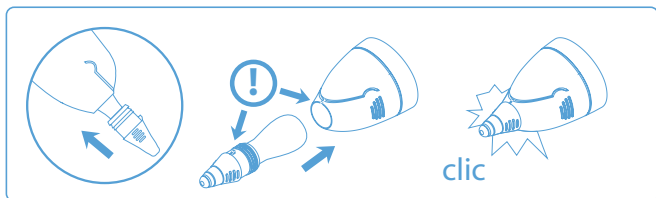

3) Position the disposable nozzle just at the opening of the nostril. To improve the efficiency of Cleanoz®, it is recommended to irrigate nasal passages with saline solution before aspirating secretions. To optimize the nose-cleaning, it is recommended to close the other nostril with your index finger while using Cleanoz® on the other nostril.

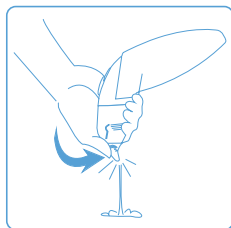


4) Switching on: Press the switch. (The secretions are collected in the single use secretion bag).

5) Switching off: The suction stops automatically a few seconds later when the bag is full.

6) Empty the reservoir (away from the baby) by pushing the nozzle tip with the thumb, before clearing the other nostril.

7) Empty, remove and discard the single use reservoir nozzle piece after use.

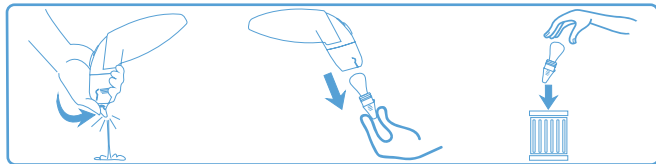

CLEANING AFTER USE



1) The nozzle piece is single use only; discard it after use.

For economical reasons, you might want to reuse the disposable reservoir nozzle piece a few hours after an earlier use.

THIS IS NOT RECOMMENDED, but if you must you should rinse the reservoir nozzle with clean water immediately after use (aspirate clear water for 3 seconds with Cleanoz® then empty the reservoir) and clean the cone shaped tip with an antiseptic cloth.

Never use the same reservoir nozzle piece on different babies!

Never reuse the same nozzle tip overnight!

2) Cleaning the main unit: Make sure that the battery compartment is shut before cleaning the main unit. Do not clean the nose cleaner with solvents or abrasive products.

3) Do not immerse the nose cleaner. Do not aspirate any other product, but water, with Cleanoz®.

4) To help protect our environment, we encourage you to recycle at approved recycling locations. Do not dispose batteries in normal household waste.



STORAGE

1) Store Cleanoz® in a dry place at room temperature. Cleanoz® and the Disposable reservoir Nozzle pieces must be stored in appropriate conditions: temp between 41° and 113°F (5° and 45°C).

2) If not used for several months, it is recommended to remove the batteries.

PRECAUTIONS

1) Cleanoz® is not a toy. Cleanoz® must be used by adults after reading this manual.

2) Keep Cleanoz® and Nozzle pieces out of children's reach.

3) Cleanoz® is designed for household use, only to aspirate baby's nasal secretions.

4) Cleanoz® is designed for an intermittent use. (15 seconds/minute max). Using the appliance without a nozzle piece or keeping the switch on when the reservoir bag is full will use batteries prematurely.

5) Before using Cleanoz®, make sure that all parts are correctly assembled.

6) Respect the expiration date of the disposable reservoir nozzle pieces.

7) Do not attempt to disassemble the nose cleaner unit or the nozzle pieces.

8) Do not strike, bend or drop the unit.



9) Use only Ubimed® single use reservoir nozzle pieces for Cleanoz® MB002. They are designed to be used only with Cleanoz® MB002

10) The Reservoir Nozzle piece is not sterile.

11) Do not use or store the unit or the disposable nozzle pieces under the following conditions:

- *Hot locations such as the inside of a car parked outside in hot weather.

- *Under direct sunlight.

- * Humid or corrosive locations.

- *Keep away from hot surfaces, inflammable substances or explosive gas

12) The nose cleaner should not be used in any conditions other than those described above .

WARRANTIES AND SERVICES

A LIMITED WARRANTY FOR THE TIME STATED ON THE PACKAGE IS PROVIDED FROM THE DATE OF PURCHASE TO THE CONSUMER WHO BOUGHT THIS PRODUCT.

WARRANTY RESTRICTIONS: WHAT IS NOT COVERED UNDER WARRANTY:

DISPOSABLE NOSE PIECES, BATTERIES. DAMAGE CAUSED BY MISUSE, ABUSE, NEGLIGENCE OR ALTERATIONS. NORMAL WEAR AND TEAR, INCLUDING CHIPS, SCRATCHES, ABRASION, DISCOLORATION OR FADING.

IMPLIED WARRANTIES:

ANY IMPLIED WARRANTIES INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH IN THIS LIMITED WARRANTY. SOME STATES DO NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE LIMITATION MAY NOT APPLY TO YOU.

LIMITATION OF REMEDIES:

IN NO EVENT SHALL UBIMED®, Inc. OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON MISUSE, BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE. LOSS OF PROFIT. LOSS OF USE; THE CLAIMS OF THIRD PARTIES AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OF LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

PURCHASE DATE: SERIAL NUMBER:

RETAILER STAMP



FCC part 15 compliant, IP 40



Components **must be recycled**



Use 3 AAA LR03 Alkaline Batteries 4,5V DC

Protected by US and International Patents Pending. Design US Patent 542 915 .
UBIMED® and CLEANOZ® are registered trademarks of Ubimed®.

UBIMED Inc, Los Angeles, California 90035 US
www.ubimed.com call center: 1-866-5-UBIMED
2007Ubimed® Inc. All rights reserved.